

(51) Int. Cl.⁶: B01J 19/00

(22) Date of filing: 21.03.1997

- **Duffield, Mark J.**
Cambridge, Massachusetts 02139-4809 (US)
- **Hogan, Eleanor M.**
Ann Arbor, Michigan 48105 (US)
- **Nickell, David Glenn**
Ann Arbor, Michigan 48103 (US)

(74) Representative: Mansmann, Ivo et al
Warner-Lambert Company,
Legal Division,
Patent Department,
c/o Gödecke AG,
Mooswaldallee 1
79090 Freiburg (DE)

(72) Inventors:
• **Dewitt, Shella Helen Hobbs**
Canton, Michigan 49236 (US)

(57) A system for developing and performing multiple simultaneous synthesis with an array of reaction vessels includes an information management system, an equipment controller interface and responsive equipment such as a liquid handling robot. Equipment may also be supplemented by other reaction peripherals and pick and place robotic equipment controlled by an event scheduler. The preferred information management system includes a design module accessing reagent, corporate chemical and reaction databases, as well as a core module that includes a reaction workbook. The

reaction workbook identifies the products to be produced, preferably by the product string of chemical building blocks used to form the product, the solid or solution parameters involved to generate the products, and the selection of reagents required to carry out the reactions. In addition, the reaction workbook organizes the reaction date to guide the equipment controller, preferably for movement relative to an array of multiple reaction vessels in a reactor apparatus.

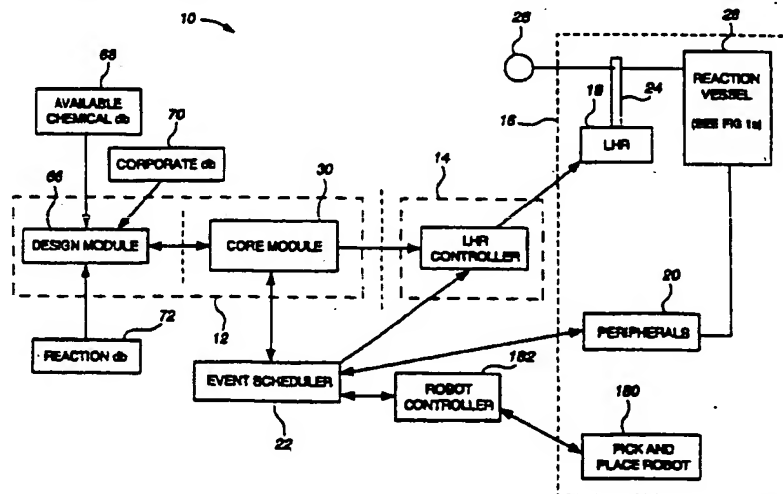


FIG-1

Description

TECHNICAL FIELD

The present invention relates generally to a system for use in multiple simultaneous synthesis of compounds, particularly an information management system for automating the preparation of chemical compounds utilizing chemical synthesis technology.

BACKGROUND ART

There are various methods and apparatuses known today for use in the multiple simultaneous synthesis of compounds, including organic compounds. For example, it is known in the art that peptides and oligonucleotides may be multiply and simultaneously synthesized utilizing various methods and procedures. It is also known that general organic compounds can be synthesized, although the procedures and apparatuses utilized are often quite different and more versatile than with synthesis of peptides and oligonucleotides. In this regard, the synthesis of general organic compounds often requires such varied conditions as utilization of an inert atmosphere, heating, cooling, agitation and/or an environment to facilitate reflux. Preferred methods and apparatuses for the multiple simultaneous synthesis of organic and other chemical compounds are set forth in U.S. Patent No. 5,324,483.

Chemical synthesis to generate compounds having possible suitability as stable, orally active drug candidates is a highly active field today. There are several methods and procedures which are currently being used to synthesize and/or screen as many chemical compounds as possible to search for new drug candidates. In order to determine the effectiveness of certain compounds in this regard, it is an advantage to be able to produce and screen thousands of chemical compounds as rapidly and accurately as possible. It is also advantageous to perform the synthesis and other processes with minimum manual interaction, and thus with as much automation as possible.

Multiple simultaneous synthesis (also called "parallel synthesis") is typically carried out in an array format which is compatible with standard techniques of organic synthesis. The final compounds are produced individually in soluble form typically on solid supports (for example, resins), although solution techniques can also be utilized. The growing compounds are readily separable from the by-products and reagents. The chemistries utilized encompass a wide variety of organic reactions.

In accordance with U.S. Patent No. 5,324,483, solid supports are contained within a plurality of gas dispersion tubes and the synthesis techniques enable the multiple, simultaneous synthesis of, for example, 8, 40, 100 or more reactions at one time.

The steps necessary to perform a synthesis generally comprise the development of a synthetic route that will be feasible on a solid support or in solution, the ver-

ification of resin-based or solution-based synthesis, and the execution of the multiple, simultaneous synthesis within an array format. The method generally involves the sequential coupling of chemical building blocks to form resin-bound or solution-based intermediates until the final compound at each array location is constructed. In addition, coupling agents or reagents are often required to help form the bond between the solid support and the building blocks.

The speed and accuracy of many of the synthesis procedures can be improved through use of a liquid handling robot ("LHR"), such as a robotic sample processor. Much of the information used to direct and perform the synthesis procedures and operate the LHR, however, is generated manually or through use of certain computerized programs. This is still labor intensive and requires significant time and effort.

In particular, the multiple reaction steps and physical processes must be determined and the proper constituents introduced and combined for each position in the array. Moreover, control systems for robots that could transport or mix the reaction compounds at each of the array position locations are subject to logic commands substantially different than the commands which are used by the product technician to create the products to be developed and tested at each position of the array. As a result, known systems have not succeeded in automating the production of testable products in an array without substantial cross training between the technical disciplines of product development and robot control to reduce product creation time.

SUMMARY OF THE INVENTION

The present invention provides a system for automating multiple simultaneous synthesis of a plurality of compounds, preferably to determine their potential as possible new drug candidates. The present invention is not limited for use only for drug applications, however, and it can be used in any industry which uses chemical synthesis and reactions to develop, test and/or screen new chemical products. For example, the present invention can be used to develop and screen new agricultural agents, dye products, diagnostic reagents, food-related agents, and the like. Also, as to the drug field, the new drug candidates can be for either human or veterinary use. An additional alternative use of the invention is for the optimization of chemical reactions, for example, reaction yields or reaction times, on solid supports or in solution.

The compound development takes place in an array, preferably provided by a reactor apparatus as disclosed in U.S. Patent No. 5,324,483 and incorporated herein by reference, although other known or equivalent reactor apparatuses can be utilized. The reactor is driven automatically by chemical information software integrated with interfacing robot control transformation software to external equipment, for example, robots and reaction implementation modules. The invention expe-

rites synthesis of products with less complications and in a quicker, easier and more efficient manner than with previously known synthesis systems.

Preferably, the system utilizes a method and apparatus for automated synthesis which comprises a plurality of modular components that automate the sequential coupling of building blocks to form intermediates until a production process is performed to produce a compound at each reaction vessel in the array of a reactor apparatus. The coupling of building blocks, the addition of coupling agents or reagents, and the manipulation of the reactions (such as by heating, cooling, agitation, etc.) are implemented automatically by equipment such as liquid handling robots or other robotic apparatus, responsive to the sequential coupling reactions chosen by a trained practitioner. The final compounds are cleaved from the spent resins to yield products for biological assay testing in substantially less time and with more efficiency than would otherwise be required to generate a family of products to be tested.

As a result, the present invention provides a system for providing multiple reaction vessels for the reactions to be performed in the production of test samples for the rapid screening of multiple chemical entities in high through-put biological assays. The present invention also provides software for orienting the practitioner's focus on chemical transformations in automated parallel synthetic reactions rather than on spatial position control of the robotic apparatus associated with a product generating station. As a result, the present invention provides substantially more comprehensive testing and precise creation of desired products than previously known methods including those employing previously known multiple simultaneous synthesis apparatus. The present invention also substantially reduces the amount of documentation and calculations previously needed or utilized by the synthesis practitioner.

In the preferred embodiment, the method includes the steps of creating an array formatted to coincide with a plurality of gas dispersion tubes and reaction wells in a holder block of a reactor apparatus. The synthesis practitioner selects the reagents, the selection becoming interfaced with a source location database with identification of particular reagent solutions and building blocks. By selecting the product structure desired, the resin weight to be provided at each reaction vessel and the product string describing the product, the present invention provides the selection and automatic distribution of reagents and processing steps to be sequentially provided at each reaction vessel. Accordingly, by concentrating on the chemical transformations to be accomplished at each reaction vessel, the present invention automates operation of the external modules without imposing numerical control or robot displacement characteristics on the commands given by the practitioner.

BRIEF DESCRIPTION OF THE DRAWINGS

The present invention will be more clearly understood by reference to the following detailed description of the preferred embodiment(s) when read in conjunction with the accompanying drawings and appended claims:

FIGURE 1 is a schematic illustration of a multiple simultaneous synthesis system and a multiple parallel synthesis system in accordance with the present invention;

FIGURE 1A depicts a reactor apparatus for use with the present invention;

FIGURE 2 is a flow chart outlining operation of the core module shown in the system of Figure 1;

FIGURE 3 is a flow chart showing process steps in the initialization of a new array in the core module of Figure 2;

FIGURE 4 is a flow chart of process steps for entering reagents to create products automatically with the core module of Figure 2;

FIGURE 5 is a flow chart of process steps for entering products in the core module of Figure 2;

FIGURE 6 is a flow chart of process steps for entering resin parameters in the core module of Figure 2;

FIGURE 7 is a flow chart of process steps for defining reaction parameters in the core module of Figure 2;

FIGURE 8 is a representative screen diagram including an overlay for dialogue with the core module of Figure 2;

FIGURE 9 is a representative screen diagram for a worksheet used in dialogue with the core module of Figure 2;

FIGURES 10a and 10b depict a worksheet similar to Figure 9 but showing additional interface capability for dialogue with a core module of Figure 2;

FIGURE 11 is a representative screen diagram similar to Figures 9 and 10 but showing additional interactive capability for dialogue with the core module of Figure 2; and

FIGURE 12 is a representative screen diagram of a worksheet displaying an instantaneous condition in an array corresponding to a physical reaction of the reactor apparatus for a system shown in Figure 1.

BEST MODE(S) FOR PRACTICING THE INVENTION

Referring first to Figure 1, the automated multiple simultaneous synthesis system 10 in accordance with the present invention comprises an information management system 12, an interface 14 and physically interacting accessories or equipment 16. In this regard, as used herein, the terms "multiple simultaneous synthesis" and "parallel synthesis" are used interchangeably and refer to the same general system and procedures. On the other hand, the term "multiple parallel synthesis" refers to a system wherein two or more

multiple simultaneous synthesis procedures or parallel synthesis procedures are part of an overall system; this is discussed in more detail below with reference to Figure 1 and the use of a "pick and place" robot.

Components in each of the units 12, 14 and 16 in the system 10 are modular to allow individual components to be substituted as needed. For example, in the preferred embodiment, the interactive equipment 16 comprises a liquid handling robot (LHR) 18, for example, a Cyberlab Model C200, although other comparable and equivalent LHR systems could be utilized, such as a Tecan robotic sample processor.

Peripherals 20 provide the capability for more complex processing or coordinated concurrent parallel synthetic reactions operated through an event scheduler 22 to improve utilization of each component in the accessories module 16. Module 16 is described herein for the sake of brevity in terms of the liquid handling robot 18 and its control of the manipulator 24, for example, interfacing between a chemical supply 26 and a reactor or synthesis apparatus 28 used in performing multiple simultaneous synthesis. It will be understood, however, that the exemplary discussion is not limiting and does not restrict the operations that can be performed by the accessories module 16.

Interface module 14 provides a control logic transformation from a core module 30 in the system 12 to the equipment and accessories. For example, the information compiled in the core module 30 may be transformed to computerized numerical control displacement paths of the manipulator 24 (such as the arms of the LHR) that define movements of fluid supply needles on the manipulator 24 that correspond to the positions of the reaction vessels of an array 56 (see Figure 1A) in a reactor or synthesis apparatus 28.

The reactor apparatus 28 preferably is of the type disclosed in U.S. Patent No. 5,324,483, which is concurrently owned by the same assignee as the present invention. The disclosure of the '483 patent is hereby incorporated by reference herein. Other reactor apparatus and synthesis systems could be used with or make use of the present invention, including those of Advanced Chemtech, Affymax, Biotech Instruments, Chiron Corp., Gibson/Abimed, Lindsey, Protein Technologies/Rainin, Rhone-Poulenc, Shimadzu, Takeda, Fuchs/PASS and Zinsser.

In accordance with the '483 patent, and as generally shown in Figure 1A, the reactor or synthesis apparatus 28 preferably used with the present invention includes a plurality of reaction tubes 32 that are each received in a reaction well 34 carried in a reservoir block 36. The lower end of each reaction tube 32 includes a filter 38 and one or more pressure release openings 40. Preferably, the reaction tubes are glass tubes with sintered glass filter bottom portions (also called "PINs") that hold resin-bound building blocks of a product. As used herein, the term "reaction vessel" refers to the physical member or structure in which the chemical synthesis takes place or is carried out for each product to

be developed in the reactor apparatus. For example, in accordance with the disclosure of the '483 patent, a "reaction vessel" includes both a PIN (reaction tube) and a reaction well.

The reactor apparatus 28 includes a holder block 42 with a plurality of apertures adapted to receive the reaction tubes 32 above the reservoir block. The apparatus 28 also includes a manifold 44 having an opening at the lower end to receive the array 56 of reaction tubes 32 while the top wall of the manifold 44 contains apertures that register with openings in the reaction tubes 32. The top 46 of the manifold 44 is covered with a gasket 48 held in place by an apertured plate 50 providing a covered access opening to each of the reaction tubes. For example, a fluid carrying needle of the manipulator 24 may penetrate the self-sealing gasket 48 and supply or remove reagent or other materials within the reaction tube 32. In this regard, the gasket may be made of rubber or a rubber-type material (e.g., silicone) which is puncturable with a needle-like object and automatically resealable following the puncture.

Ports 52 and 54 in manifold 44 provide an inlet port and an outlet port, respectively, for changing the atmosphere about the reaction tubes 32. As a result, the reactor 28 provides an array 56 of reaction wells, that can be subjected to various process steps such as agitation, heating, cooling, etc. as desired, in order to conduct multiple simultaneous synthesis as more fully discussed in U.S. Patent No. 5,324,483. Connectors, for example, clips 58, secure the reservoir block 36 to the holding block 42 about a sealing gasket 60, while clips 62 secure the holding reservoir 42 to the manifold 44 about a sealing gasket 64. Fasteners are also used to hold the apertured plate 50 to the top of the manifold. Other mechanisms for holding together the components of the reactor apparatus could also be utilized.

The information management system module 12 comprises a design module 66 interfacing with the core module 30. The design module 66 enumerates a series of compounds based on a common core molecule and a given set of substituents available from a plurality of databases. For example, an available commercial chemical database 68, such as Available Chemical Directory (ACD) by Molecular Design Limited (MDL), a corporate database 70 of chemical structures, and a reaction database 72, such as REACCS by MDL, or Cross-Fire by Beilstein, can be used to define a library of compounds in the design module 66. A corporate database, for example, is data or information accumulated or developed by the company or entity which is utilizing the invention. These compound reactions and reagents are then exported using a common information exchange format, for example a structure data file (SD), or reaction data file (RD) to the core module 30 for use in a "reaction workbook". The chemical directory in a commercially available database, such as the Maybridge database, provides the decisions regarding which compounds are to be made based on several criteria such as availability and cost of the reagents. The

design module 66 integrates the information contained in the databases including the particular compounds already available in the corporate repository or database.

As indicated, the reaction workbook of the core module 30 interfaces with the reaction based databases such as MDL's REACCS and Beilstein's Cross-Fire in conjunction with commercially available chemical databases to provide the basis for initializing and performing chemical reactions by the automated synthesis system 10. Data about conditions and generic reagents are downloaded from these databases to the core module 30 via the design module 66. Several commercially available products embody some of the functions of a design module and are incorporated as desired. For example, MDL's Project Library, Tripos' Legion and Selector, the DiverSelector software and the CombinD-BMaker may be used as the design modules.

The core module 30 simplifies the organization and coordination of multiple simultaneous synthesis in the array 56 by transforming reaction data to define destination data and maps that correspond to process steps performed in the array 56. The destination data is transported to a robot controller by a defined interface for data transfer, for example, ASCII or SD files, so that the robot is driven by controller data structured according to specific process steps in the array. While the core module 30 may generate a reaction workbook with the transported data in response to downloaded product and reaction data from a design module, the core module may also create a reaction workbook by independent user interaction with a graphical user interface.

The graphical user interface simplifies the practitioner's interaction with the equipment 16 and the available database information. The system user inputs data of the type previously written in a laboratory notebook into the reaction workbook of the core module 30. This data includes lists of the products, the amounts, the number of mols of constituents, the number of equivalents for each reagent, the molecular weights for each reagent, and a description of the procedure to be performed. This data generates control instructions for the robot controller software. At the end of the synthesis, the results of the reaction are recorded, including the percent of yield, the weight of material obtained, and the purity of the material.

In order to drive an automated chemical synthesizer, the reaction workbook formalizes the identification of the scheme for synthetic transformation, the identification of the constituents involved, and the identification of the locations at which each synthesis sample is being transformed. The reaction workbook automatically assigns a destination in the array to each product entered. The key to linking the products with the reagents and the positions in an array is called the "product string". The product string is a concatenation of the reagents used to construct the desired compound and placed in the order that the reagents are used in the reaction sequence and chemical transformation.

The reagents may only be a part of the target molecule temporarily, for example, the resin used in a solid phase reaction. For example, a product string AA,BB:CC:DD would describe the molecule AABBC-CDD formed from a three-step reaction. The first step introduces the building blocks AA and BB. As another example, in a four-step solid phase reaction sequence in which a reagent is attached to the resin in the first step and the target molecule is cleaved in the last step, the product string will be AA:BB:CC: None.

In addition, the reaction workbook programming permits the grouping of reagents to improve the efficiency of reagent handling and sorting and reduces time in a logical manner. For example, if two building blocks are used in the first step of synthesis in an array, the array can be logically divided in half with all the reaction vessels in one half receiving one building block while manipulator 24 provides a second building block to the other half. In addition, the array grouping may be a pattern identified by a destination map generated by the reaction workbook.

The reaction workbook can also accommodate different chemical reactions in a single synthesizer unit by linking reagents and solvents to specific building blocks. For example, both acid chlorides and organic acids could be used simultaneously in the same array to introduce different building blocks by linking the necessary reagents to either the acid chloride or the acid. Also, either solid supports (typically a resin) or standard solution techniques can be utilized in the synthesis procedures in accordance with the present invention.

The reaction workbook is implemented in a spreadsheet paradigm, preferably with the chemical structures that graphically represent molecule formation and other information contained in cells. The Product Worksheet will have a place to specify a generic structure and columns for groups to be modified. Molecular weights and formulas can be automatically calculated from the structures of products. Procedures for each reaction of a multiple step sequence can be placed on a separate Reaction worksheet that includes both the representation of the reaction and the chemical structures. Based upon the amount of starting material input, the Reaction spreadsheets will calculate the amount of each reagent needed for each step. For solid phase reactions, the average weight of resin in each reaction vessel along with the resin loading estimation of reactive sites, measured in equivalents per weight, will be used to calculate the number of mols of starting material. Alternatively, the program allows the input of individual weights to be assigned to each reaction vessel. Finally, the program will check to see that the total volume of reagents and solvents added at an individual PIN does not exceed the allowable volume of the reaction well of the array 56.

Referring to Figure 2, the general format 74 of the reaction workbook in the core module 30 is shown utilizing input from a design module 66. After a choice of operating levels 75 or 76, spreadsheets 77 or dialogue boxes 78 are used to identify an existing array 80 or to

create a new array 82. The user may choose to enter information into the reaction workbook or transport data to the workbook from the design module. When the workbook is open, the user may enter from a menu Reagents as shown at 84, Products as shown at 86, Resin Parameters as shown at 88, and Reaction Parameters as shown at 90. Based on this data, the Export Synthesis Procedure as shown at 92 is generated.

Preferably, the program functions on the menu are accessible from a graphical user interface employed with a computer terminal work station. For example, in the preferred embodiment, Microsoft Excel provides templates developed for accessing functions of the system according to the present invention. The functions of Create New Array 82 (as shown in Figure 3), Enter Reagents 84 (as shown in Figure 4), Enter Products 86 (as shown in Figure 5), Enter Resin Parameters 88 (as shown in Figure 6), and Enter Reaction Parameters 90 (as shown in Figure 7) can be accessed directly from the main menu, for example, by the use of graphic buttons on a computer screen. In addition, steps performed in each of these functions can be accessed through spreadsheet templates identified at the screen portions shown as tabs along the bottom edge of the screen. The buttons or functions can be selected or accessed by use of a touch screen, by the clicking of a mouse, or by the use of a keyboard, depending on the type of computer equipment utilized.

With reference to Figure 8, graphic buttons 94, 96, and 98 are provided on the computer spreadsheet diagram 93 in order to access the Create New Array, Enter Products and Reaction Parameter menus, respectively. The main menu 93 also contains graphic buttons for Create Reagent List and Resin Parameters, but they are covered by the dialogue box shown. Also tabs 100, 102, 104, 106, 108, 110 and 112 on the screen diagram can be used to access worksheets relative to Products, Reagents, Resin and the like.

Upon actuation of Create New Array graphic button 94, a dialogue box requests the user to enter the number of reactions to be involved in the creation of a product, and displays autographic buttons to select the size of the array. In the preferred embodiment, the dialogue box is set up with a selection between an eight PIN and a forty PIN array, although it is to be understood that arrays other than PIN arrays and with different numbers of reaction vessels can be utilized depending upon the products and the reactor apparatus to be utilized. A title for the array is entered and the work product is saved after the worksheets for the reaction workbook are initialized, and the array can be developed using the autographic button functions for worksheets.

Referring now to Figure 4, the Enter Reagents 84 function of the preferred embodiment includes entering the reagent name 200, reagent code 202, molecular weight of the reagent 204, density or concentration of the reagent 206, preferably molarity, and the type of reagent used 208. A spreadsheet 105 showing this data is

shown in Figure 9. In a preferred embodiment, the reagent type is selected from: (a) building blocks that build upon the resin carried in the reaction vessel, (b) solvents that carry or affect but do not join the building blocks, (c) reagents that participate in a reaction but do not end up in the target molecule, (d) resin-bound starting material, and (e) resin-bound reagents. The ease with which reagent selections can be made is shown in Figure 8 in which a dialogue box 114 queries the user for selection of the appropriate constituent. In addition, the box 114 provides a listing window 116 in which a cumulative list of the selected reagent names for an array are added or deleted. The window list corresponds with the reagent data stored in the reagent worksheet 105, depicted by the tab 104 in Figure 8. The source location information is automatically assigned by the program.

The Enter Products function 86 of the preferred embodiment shown in Figure 5 similarly enables the user to select a product name, the molecular formula, and the building blocks for each transformation. In addition, the array location corresponding to the PIN at which the product is to be made is automatically assigned. Also, the product string is entered, and the molecular weight is automatically calculated and stored in the Product worksheet accessible at the lower edge tab 102 in Figure 8. In addition, the product strings are parsed to identify the discrete building blocks and stored. The products may then be sorted by recognition of the unique building blocks required at each array location to build a product at each reaction vessel of the reactor apparatus. Other portions of the worksheet, for example, the chemical structure, can be copied and pasted from another source.

As best shown in Figures 10a and 10b, the Product worksheets 103a, 103b of the preferred embodiment enables entry and listing of the product name to be defined at each PIN. In the eight PIN array shown, the worksheet also provides the listing of reagents R1 and R2 at 118 and 120 and the product strings 122 identifying the process steps that are to be performed, for example, the three process steps identified, when Creating New Array 94 has been selected for initializing an array. In addition, the following data is determined: the molecular weight 124 of the product identified by the user, the molecular formula 126 of the product entered, and the theoretical product weight 128, if 100% of the reaction occurs. In addition, a separate listing of each of the building blocks used in the first 130, second 132 and third 134 chemical reactions to occur in each step of the process creating the product is recorded (this is shown in Figure 10b which depicts spreadsheet 103b which is the second half of spreadsheet 103a). The parsing of product strings, sorting of product strings and the creation of product strings can be accessed by graphic buttons 140, 142 and 144 on the Product spreadsheet 103a, as shown in Figure 10a.

The product strings must be entered into the system to relate a specific product to the building blocks

used to create it. The codes for each building block are concatenated into a string of characters in the order they would appear in the synthetic sequence. Individual building blocks are separated by colons. As explained above, if building blocks "A", "B", "C", for example, are to be incorporated into a final product in a three-step synthetic transformation, the product string would be A:B:C. If no building block is utilized in the last step (e.g., the removal of a protecting group or the cleavage of a product from the supporting resin) the word "None" is substituted (A:B:None). If two building blocks are introduced in a single step, the individual components are separated by a comma (e.g., A,B:C:D).

By selecting the "Create Product Strings" button 144 on the "Products" spreadsheet 103a, another dialogue box is displayed. The user selects the "Product Name" on the dropdown list to display the names of the products previously entered in Enter Products 86. The product string created in this dialogue box will be entered into the appropriate cell in the "Product String" column. This dialogue box is the manual alternative to product string construction. The product string is then created by sequentially selecting the building blocks shown in a scrolling list box in the dialogue box. Product strings can also be transferred from a design module to the Reaction Workbook.

Selecting the next building block concatenates it to the product string. This interactive process is continued until all the building blocks have been selected for the product. Selecting a "Save" button in the dialogue box, enters the product string into the "Products" spreadsheet. The next product name is then selected and the process is repeated until product strings have been created for all of the products. The "Rxn1", "Rxn2", and "Rxn3" columns will be populated in the next step when the individual building blocks are parsed from the product strings.

Selecting the "Parse Product Strings" button 140 on the "Products" spreadsheet 103a creates the data shown in spreadsheet 103b of Figure 10b. The program associates each building block with a reaction step (e.g., "Rxn1") and enters it in the appropriate cell. Thus the building block "fGly" would be used in the first step of the synthesis for the product specified in row 4 of column 1.

The Reaction Workbook also contains a program which will order the reagents in an array. Selecting the "Sort Product Strings" button 142 on the "Product" spreadsheet 103a will organize building blocks by a series of sorts in the parsed building block columns. The algorithm orders the building blocks based on the number of unique building blocks in each reaction and their name. This process is useful for optimizing the placements of the reagents by a liquid handling robot or simplifying the order during manual addition. It is also possible to integrate the functions in parsing and sorting product strings into the "Product Information" dialogue box.

A multiple simultaneous synthesis can be origi-

nated by placement of a resin in each PIN in the array 56. Accordingly, the process for Entering Resin Parameters 88 is shown in Figure 6. When entering data about a solid phase organic synthesis, as shown in this example, the end user must describe the resin's characteristics to the workbook. Clicking on the "Resin Parameters" button on the main screen will display the dialogue box. Depending on the size of the PIN array, the Resin Parameters dialogue box will show a corresponding number of fields. The user may either enter the individual weights of the resin used in each PIN or may enter an estimated average weight of the resin in each PIN. It is possible to toggle between these two options by selecting one of the two buttons in an "Average Resin Weight" box. Selecting the "Enter Individual Pin Weights" button activates the data fields for data entry. When an average value has been entered, the individual fields are grayed out indicating that data may not be entered into them, while the values displayed are copied from the "Avg. Resin Wt." field. Data about the loading factor of the resin used is also entered. Selecting the "OK" button enters this data into the Resin spreadsheet. When individual weights are entered, the average weight is calculated and stored. The theoretical product yield based on that average weight is also stored in the Resin worksheet accessible at tab 106.

The Enter Reaction Parameters function 90 of the preferred embodiment is shown in Figure 7. The dialogue box 150 for selecting reaction parameters is shown displayed over the first reaction spreadsheet 109a in Figure 11. In addition, the spreadsheet at tab 108 displays a series of short actions in abbreviated form at 152 and in expanded form at 154 relative to each product string reaction, for example, the repeated addition and stirring of a particular reagent until a predetermined amount of reagent has been introduced to a PIN in the reaction vessel at a particular location. In addition to identification of the Short Action 152 and the detailed Description 154 of the short action step 152, the spreadsheet 109a also displays the Reagent employed 156, the PIN Destinations 158, the total number of PINs 160 ("PIN count") involved in the step, a total Volume 162 of solution per reaction well, and the Equivalents 164 (parts required per complete reaction) of the reagents to be added. The temperature, time and reagent weight are also entered or provided by transport from another module.

The next step in the process of entering reaction data is to provide the Reaction Workbook with the details of the actual synthetic procedure. Selecting "Reaction Parameters" from the menu results in display of dialogue box 150 shown in Figure 11. If the "Reaction Parameter" button 98 on the "Start Here" spreadsheet 100 were selected, the same dialogue box would appear but on top of that spreadsheet 93. The first step in the "Reaction Parameters" dialogue box 150 is to select the number of the reaction step into which information will be entered. This is accomplished by selecting or clicking on the "Reaction #" dropdown list 153

then selecting it from the displayed choices. The next step is to select the desired action from the "Select an Action" dropdown list 155. Other dropdown lists and data entry fields are keyed to the action selected. For example, if "Add to List" is selected, a dropdown list for selecting the building block to be added and a field to enter the number of equivalents of the building block will appear. These fields would not be displayed if, for example, "Wait for a specified time" were chosen as building blocks and equivalents are not relevant to a pause in the reaction procedure. All fields in the "Reaction Parameters" dialogue box 150 are context sensitive in this manner. The next action would be to select a building block from the "Select a Building Block" dropdown list.

In Figure 11 the result of entering three additional building blocks can be seen. The dialogue box 150 as shown is configured to add 1.0 ml of a 25% piperidine/DMF solution. The available fields in the dialogue box are different than the previous box as additional fields are pertinent to the "Add reagent or solvent" selection. The "Link Reagent" check box 157 allows the reagent which is to be added to be linked to a previous reagent. In other words, it will be added only to the same PINs as the linked reagent or building block. Figure 11 shows the entire sequence of steps for the first step of the reaction sequence.

This spreadsheet screen also provides access to an "Go to Array Template" graphic button 190. The actuation of button 190 provides an array template identification of the reagent being introduced during the reaction in accordance with a destination map 192 as shown in Figure 12. Such an array map 192 can be printed and placed on top of a reactor apparatus array to guide manual deposit of reagents.

The transfer of data from the core module 30 to guide the robot controller's preparation of commands to the liquid handling robot 18 occurs by actuating "Generate Robot Output" graphic button 194. The reaction workbook plus random access memory (RAM) transforms the chemical location and identification information into guides for controller generated movement commands at interface 14 for the LHR 18. For example, a liquid handling robot manipulator 24 will be displaced for retrieving a selected reagent from the supply 26 (Figure 1) and delivering it to each of the PINS in the array 56 as identified by the destination map 192 developed from the design module data or instructions received from the system user.

In addition to the development and insertion of ingredients that can be accomplished through the liquid handling robot 18, it will be understood that additional automation can be used in the process of multiple parallel synthesis, that is, the multiple simultaneous synthesis (or parallel synthesis) utilizing two or more reactor apparatus. For example, equipment for handling multiple reactor apparatus 28, including pick and place robots 180 such as an ADEPT system with a robotic controller 182, can be controlled and scheduled by event scheduler 22 in conjunction with operation of the

liquid handling robot 18. (This is shown in Figure 1.) In this manner, multiple parallel synthesis can be performed in multiple reactor apparatus 28 and simultaneously handled to perform different simultaneous operations in both reactor apparatus. Likewise, the function of additional peripheral equipment such as apparatus to wash the reaction wells of one reactor apparatus, can occur while developing syntheses in another reactor apparatus. The event scheduler 22 may be similarly reconfigured similar to the reaction workbook and may also include a program with graphic user interface such as Lab View by National Instruments, Austin, Texas.

The present invention has many uses for many different applications. Although the invention is described above with particular reference to its use in identifying new drug candidates, the invention can also be used, for example, to screen for veterinary agents, other pharmaceutical agents, agricultural agents, diagnostic agents, and the like. The invention has possible use in any industry which uses chemical synthesis and reactions to develop and/or screen new chemical products. As other alternatives, the invention can be used for the optimization of chemical reactions, such as reaction yields or reaction times, for statistical experimental design strategies, and for quantitative structure activity relationship strategies.

Having thus described the present invention, many modifications thereto will become apparent to those skilled in the art to which it pertains without departing from the scope and spirit of the present invention as defined in the appended claims. For example, the present invention can be employed with solution-based chemical transformation as well as solid-based chemical transformations such as the resin-based transformation of the preferred embodiment.

Claims

1. A method for conducting multiple simultaneous synthesis in a plurality of reaction vessels, said method comprising the steps of:
 - a) creating an array of chemical transformations in locations corresponding to each of said reaction vessels;
 - b) selecting products to be developed in said reaction vessels;
 - c) assigning locations for each product in said array;
 - d) selecting constituents for developing each product;
 - e) entering data of said selected products, assigned locations and selected constituents as input to a processor; and
 - f) transforming said entered data to designation data that guides logical control of an apparatus interacting with physical positions of said reaction vessels.

2. The method as set forth in Claim 1 further comprising the step of identifying a product string.
3. The method as set forth in Claim 1 further comprising the step of positioning said reaction vessels in a reactor apparatus.
4. The method as set forth in Claim 1, wherein said interacting apparatus is a liquid handling robot.
5. The method as set forth in Claim 1 further comprising the step of positioning said reaction vessels in a plurality of reactor apparatus and said method of Claim 1 is performed substantially concurrently on each of said reactor apparatus.
6. The method as set forth in Claim 1, wherein said method further comprises the step of automating data acquisition and control of peripheral equipment.
7. The method as set forth in Claim 6, wherein said peripheral equipment comprises at least a liquid handling robot.
8. The method as set forth in Claim 1, wherein said method comprises mapping and arrangement of constituent locations for at least one step in said chemical transformations.
9. The method as set forth in Claim 8, wherein said method comprises mapping an arrangement of constituent locations for each step in said chemical transformations.
10. The method as set forth in Claim 1, wherein said step of creating an array comprises inputting a predetermined number of steps to be performed in said chemical transformations.
11. The method as set forth in Claim 1, wherein said step of creating an array comprises inputting an identification number corresponding to each of said reaction vessels.
12. The method as set forth in Claim 1, wherein said step of selecting products comprises inputting resin-based parameters and said synthesis is performed on solid supports.
13. The method as set forth in Claim 1, wherein said step of selecting products comprises inputting solution-based parameters and said synthesis is performed in accordance with solution-based technology.
14. The method as set forth in Claim 1, wherein said method comprises automating operation of peripheral equipment and introducing reagents to said reaction vessels.
15. The method as set forth in Claim 1, wherein each of said reaction vessels comprises a gas dispersion tube and a reaction well.
16. A method for conducting multiple simultaneous synthesis in a plurality of reaction vessels in a reactor apparatus, each of said reaction vessels comprising a PIN and a reaction well, said method comprising the steps of:
 - a) creating an array of chemical transformations in locations corresponding to each of said reaction vessels;
 - b) selecting products to be developed in each of said reaction vessels;
 - c) assigning locations for each product in said array;
 - d) selecting constituents for developing each product;
 - e) entering data of said selected products, assigned locations and selected constituents as input to a processor; and
 - f) transforming said entered data to destination data that guides logical control of an apparatus interacting with physical positions of said reaction vessels.
17. The method as set forth in Claim 16 and further comprising identifying a product string.
18. The method as set forth in Claim 16, wherein said method comprises automating operation of peripheral equipment.
19. The method as set forth in Claim 18 further comprising introducing reagents to said reaction vessels by said peripheral equipment.
20. The method as set forth in Claim 19, wherein said reagents are grouped for introduction to a plurality of selected reaction vessels.
21. The method as set forth in Claim 16, wherein said method comprises mapping and arrangement of constituent locations for at least one step in a chemical transformation.
22. The method as set forth in Claim 21, wherein said method comprises mapping an arrangement of constituent locations for each step in a chemical transformation.
23. The method as set forth in Claim 16, wherein said step of creating an array comprises inputting a predetermined number of steps to be performed in said chemical transformations.

24. The method as set forth in Claim 16, wherein said step of creating an array comprises inputting an identification number corresponding to each of a plurality of wells. 5
25. The method as set forth in Claim 16, wherein said selecting products step comprises inputting resin-based parameters. 10
26. The method as set forth in Claim 16, wherein said selecting products step comprises inputting solution-based parameters. 15
27. The method as set forth in Claim 16 further comprising the conducting of said method for multiple parallel synthesis with a plurality of reactor apparatus, each reactor apparatus having a plurality of reaction vessels. 20
28. A system for developing a multiple simultaneous synthesis in a reactor apparatus with a plurality of reaction vessels therein, comprising: 25
- a reaction workbook for introducing multiple product development processes, including a graphical user interface for inputting location information corresponding to an array of the reaction vessels, selecting products to be developed in each reaction vessel by defining a chemical transformation, and selecting locations for each chemical transformation in said array, said workbook developing destination data for each said reaction vessel; and 30
 - a peripheral equipment control for displacing a manipulator relative to said plurality of reaction vessels as guided by said destination data. 35
29. The invention as defined in Claim 28, wherein said peripheral equipment control comprises a liquid handling robot control. 40
30. The invention as defined in Claim 28, wherein said peripheral equipment control comprises a map of constituent locations for at least one step in said chemical transformation. 45
31. The invention as defined in Claim 30, wherein said peripheral equipment control comprises a map of constituent locations for each step in said chemical transformation. 50

FIG-1

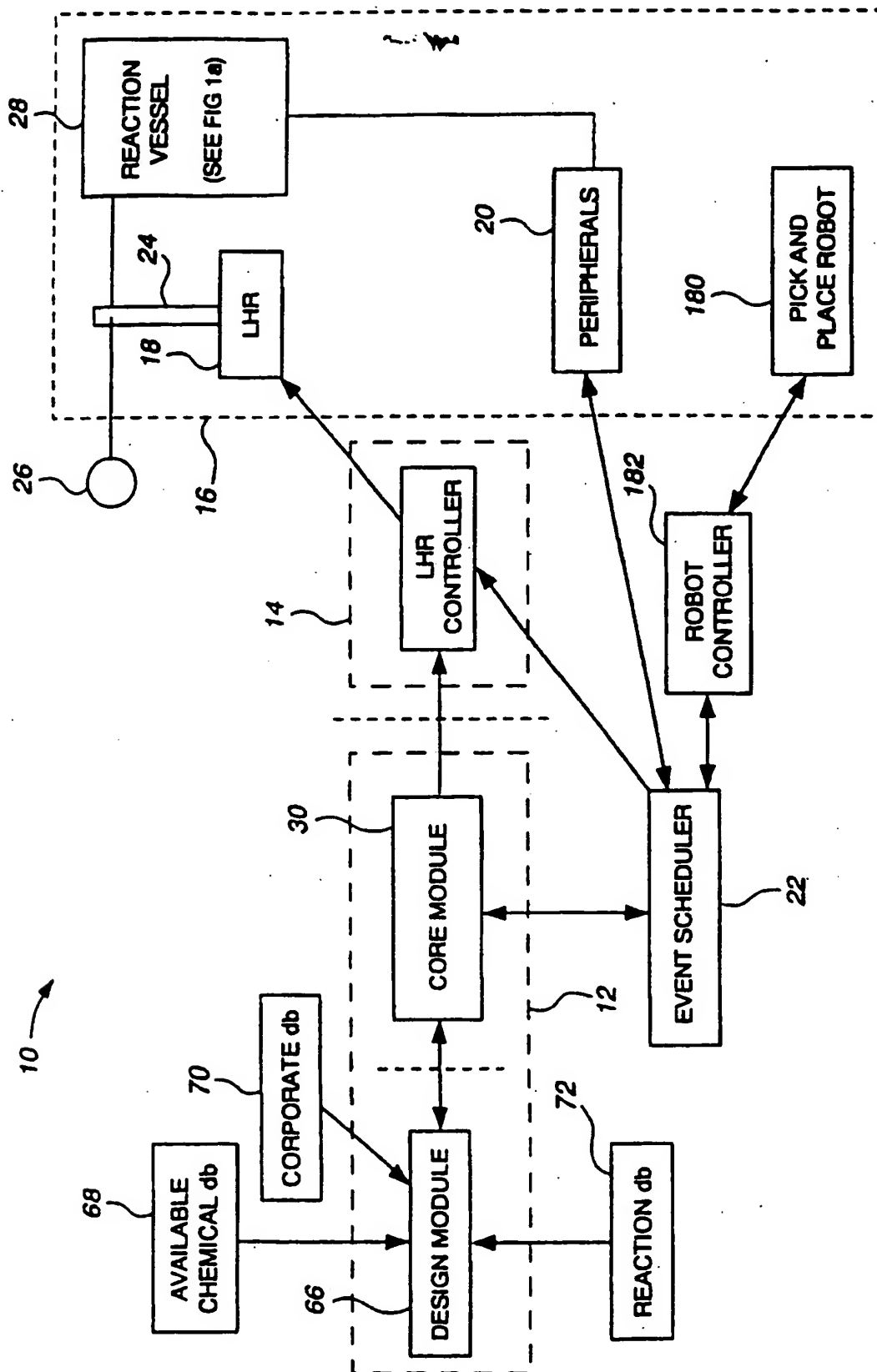


FIG-1a

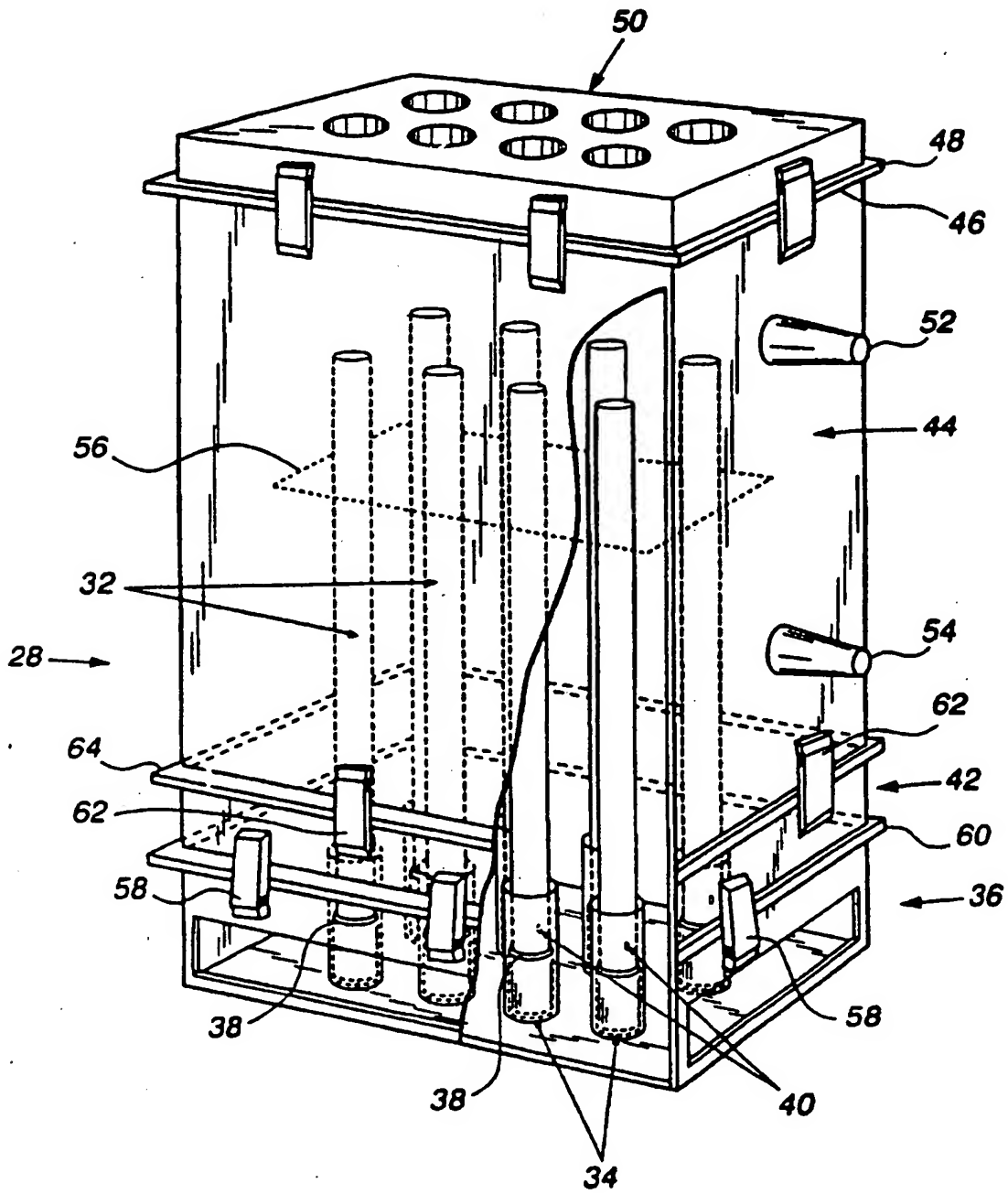


FIG-2

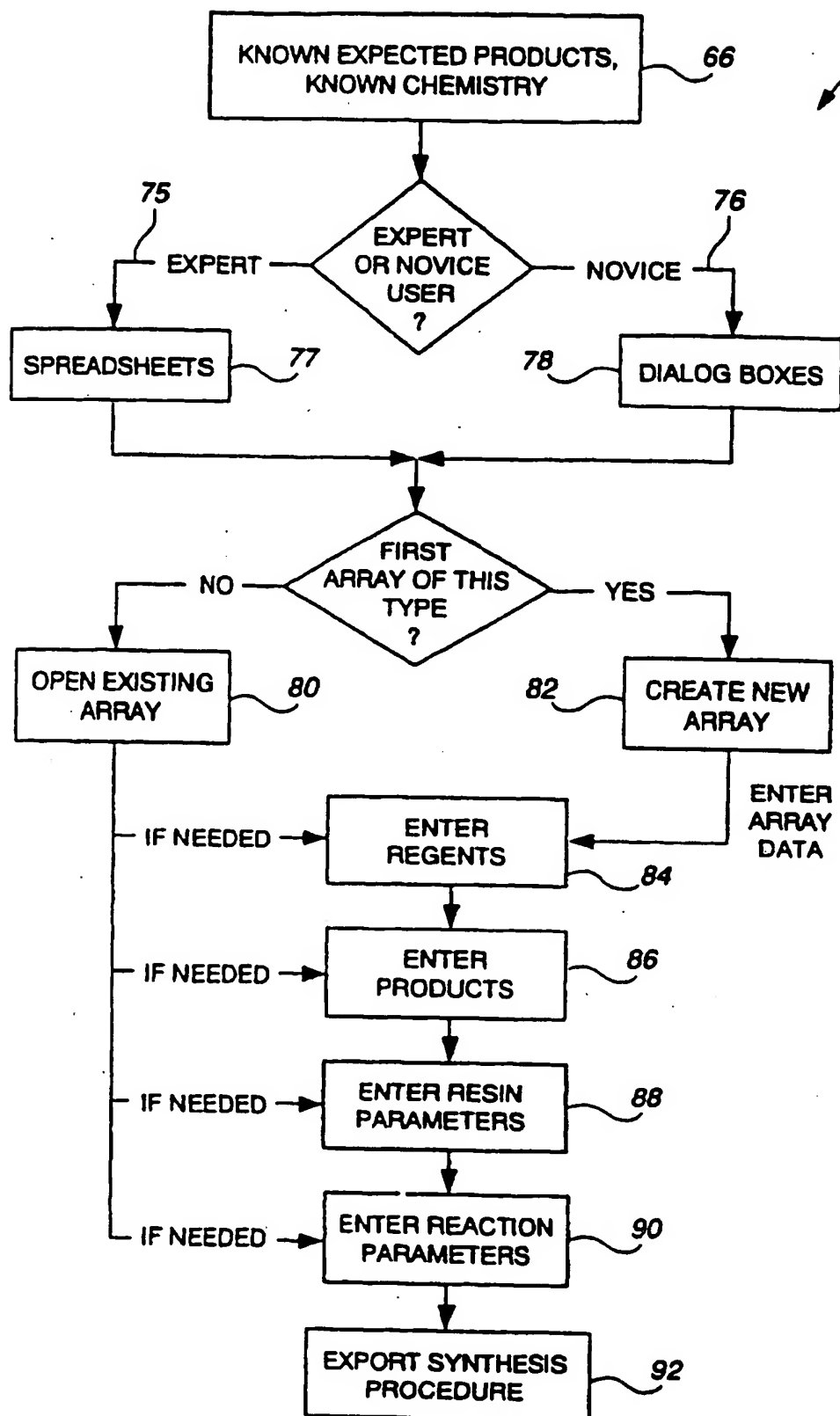


FIG-3

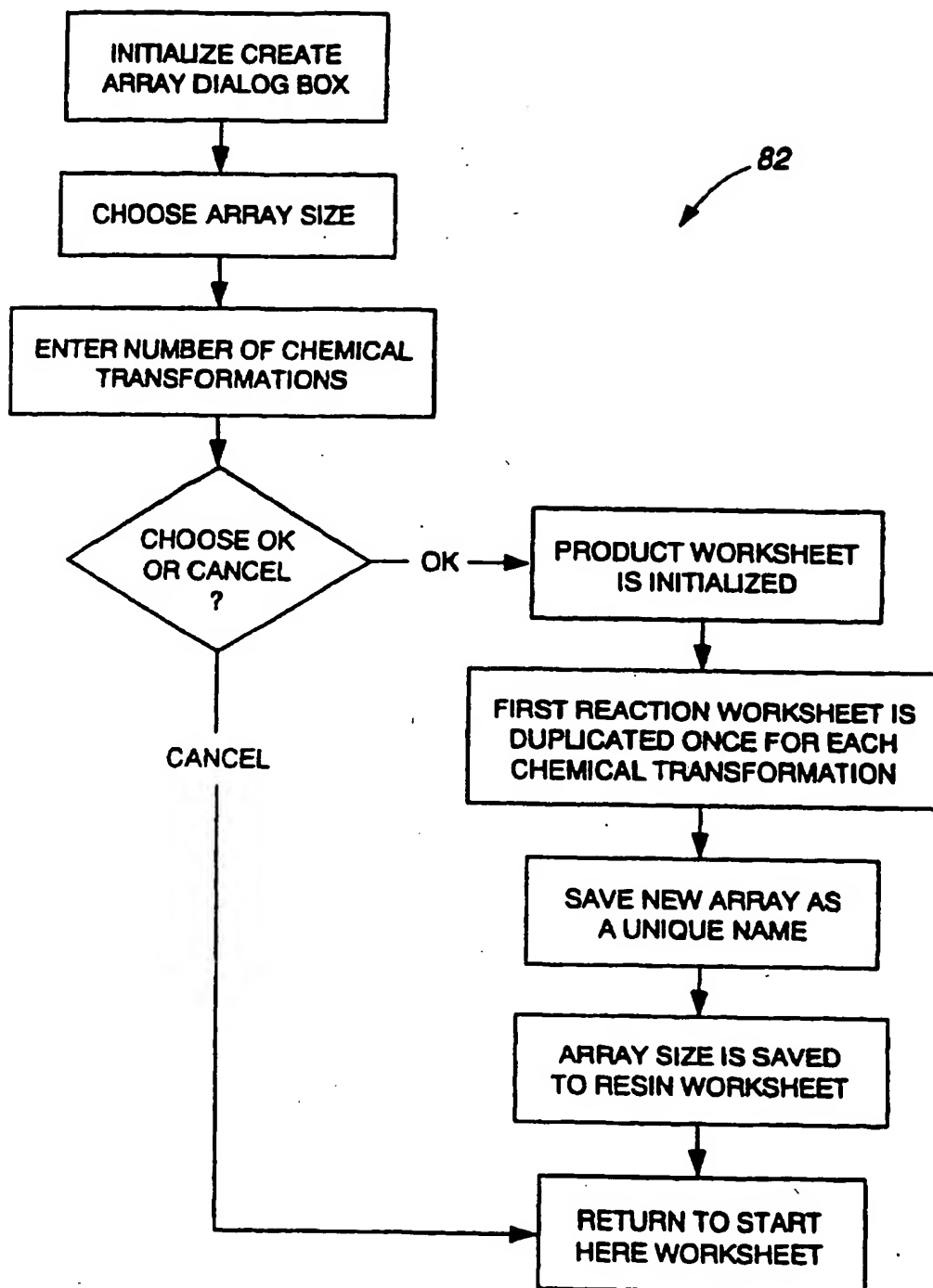


FIG-4

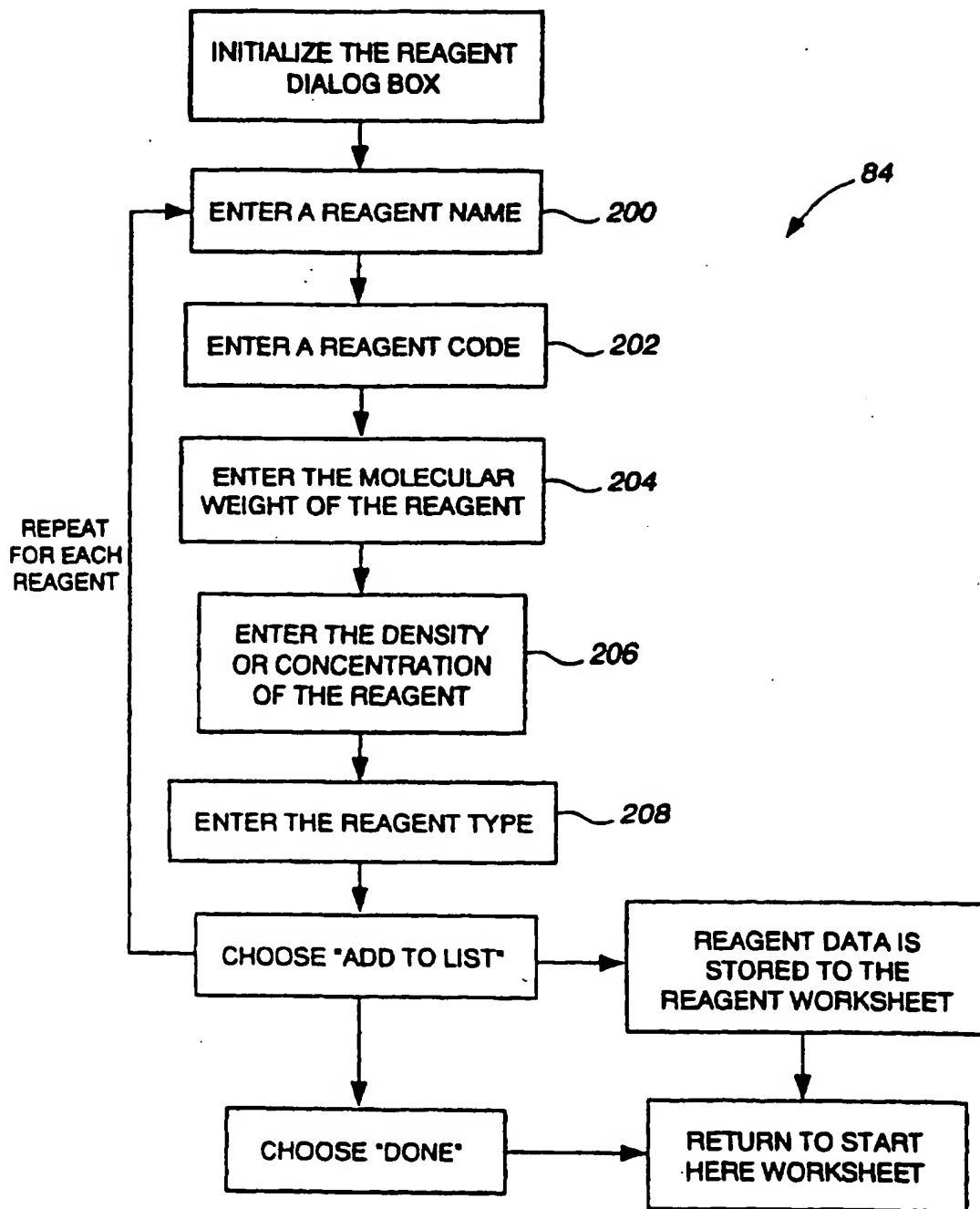


FIG-5

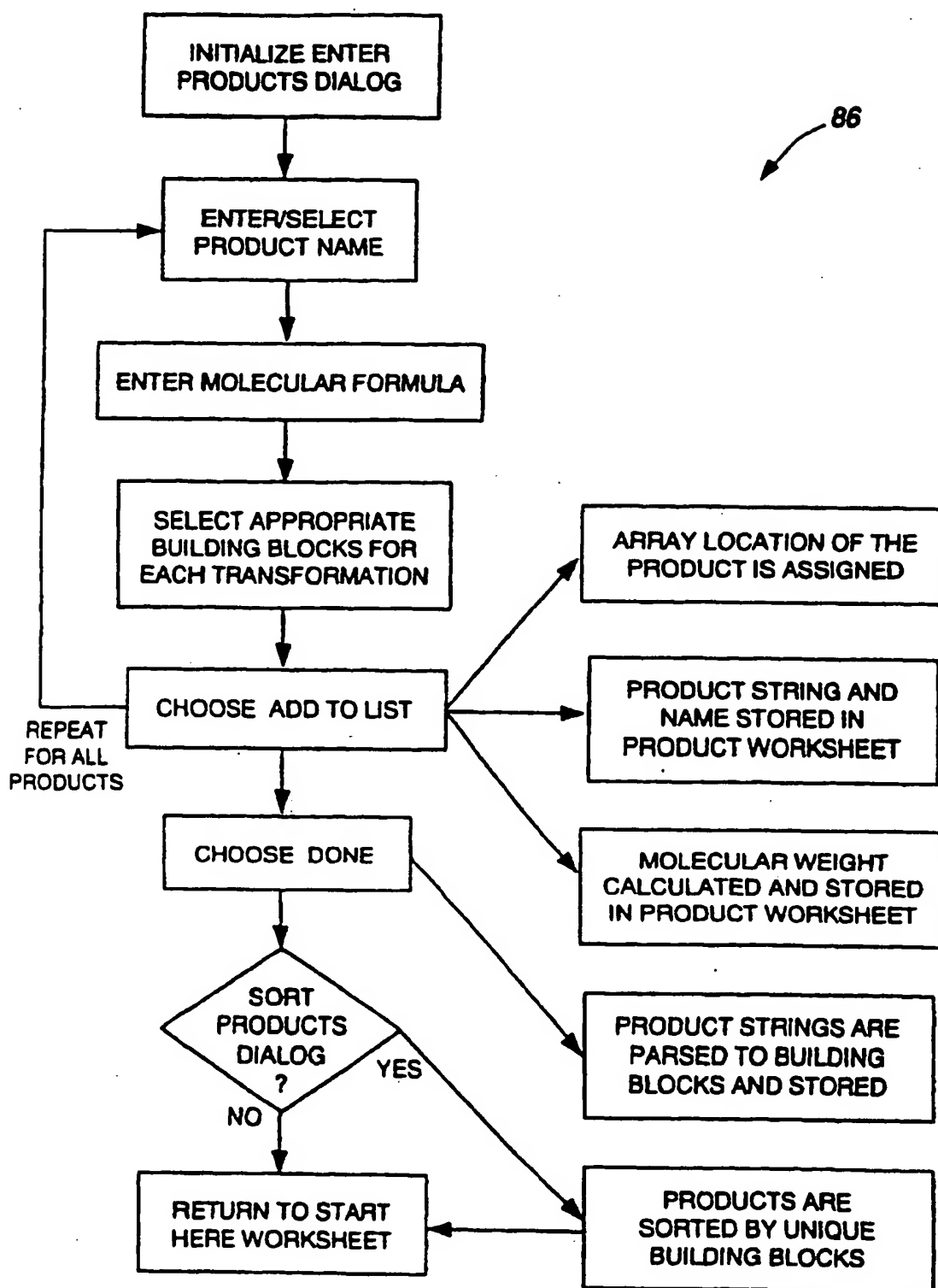


FIG-6

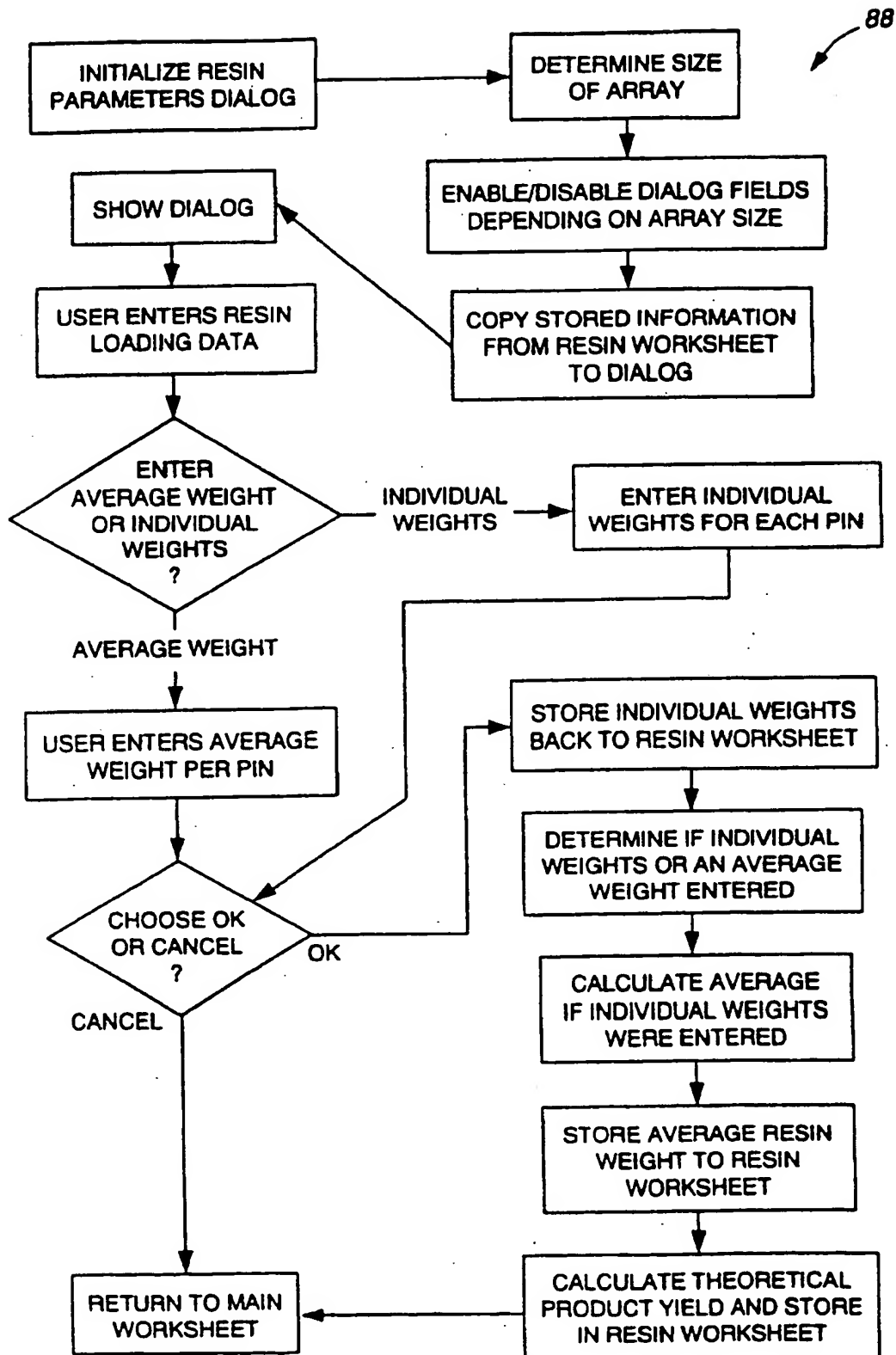


FIG-7

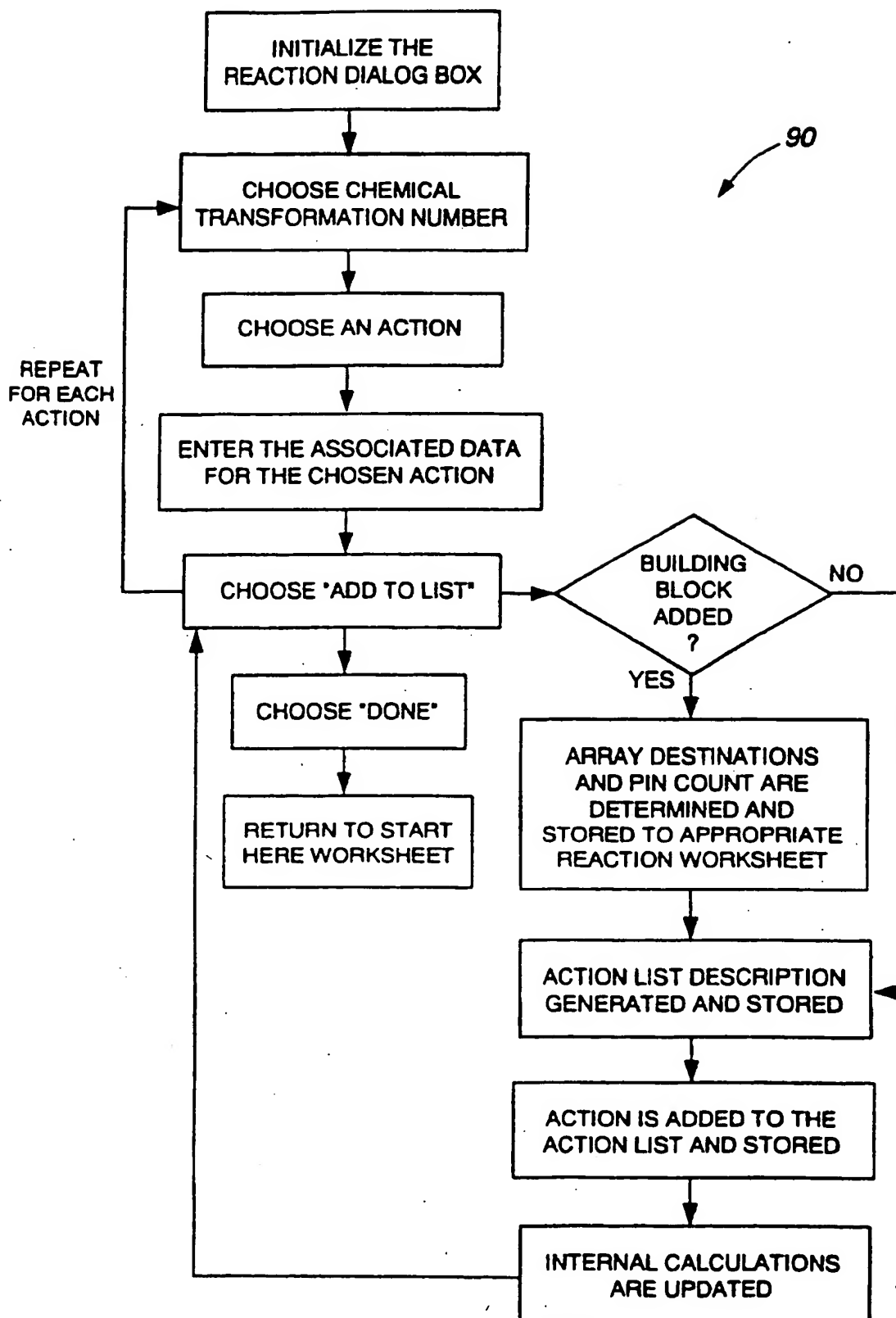


FIG-8

100 102 104 106 108 110 112

Start Here / Products / Reagents / Reactions / Reaction 1 / Reaction 2 / Reaction 3

94 Create New Array

96 Enter Products

98 Reaction Parameters

93

Reagent List

Diethyl Ether
Tetrahydrofuran
Dichloromethane

Add to List
Delete
Done

Reagent Name: Hexane
Reagent Code: HEX
Density (g/ml): 0.659
Molecular Weight: 86.18
Molarity (M):

Reagent Type:
☐ Building Block
☐ Solvent
☐ Reagent
☐ Reaction-bond Starting Material
☐ Reaction-bond Reagent

FIG-9

	A	B	C	D	E	F	G	H	I
	Reagent Name	Code	Molecular Weight	Molarity (M)	Density (g/ml)	Type	Source Location		
1	Diethyl Ether	THF	74.12		0.706	Solvent	A1		
2	Tetrahydrofuran	THF	72.11		0.889	Solvent	A2		
3	Dichloromethane	DCM	84.930		1.325	Solvent	A3		
4	Hexane	HEX	86.180		0.659	Solvent	A4		
5	N,N-Dimethylformamide	DMF	73.100		0.944	Solvent	B1		
6	Methanol	MeOH	32.040		0.791	Solvent	B2		
7	n-Butyl Acrylate	n-BuAcr	99.130		0.880	Building Block	B3		
8	P-Tolyl Isocyanate	P-TolNCO	133.150		1.056	Building Block	B4		
9	25% Piperidine/DMF	PIP/DMF				Reagent	C1		
10	Triethylamine	Et3N	101.190		0.726	Reagent	C2		
11	12N HCl	HCl		12.000		Reagent	C3		
12	Trimac Phenylamine Resin	THPa				Resin Bound Starting Material	C4		
13	Trimac Glycine Resin	THG				Resin Bound Starting Material	D1		
14	Trimac Alanine Resin	THA				Resin Bound Starting Material	D2		
15	Trimac Valine Resin	THV				Resin Bound Starting Material	D3		
16	None	None							
17									
18									
19									
20									

100	102	104	106	108	110	112
44						
45						
46						
47						
48						
49						
50						
51						

105

FIG-10A

FIN 10A

118 120 140 142 122 124 144 126 128

Hydantoin.xls

Create Product Strings

Parse Product Strings

Sort Product Strings

Molecular Weight

Theory Prod. Wt. (mg)

Molecular Formula

Product String

R1 R2

PRODUCT NAME

PM ID

3-Butyl-2,4 Imidazolidinedione

3-Butyl-2-(phenylmethyl) 2,4 Imidazolidinedione

3-Butyl-2,4-(1-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

3-Butyl 2,4-(4-methylphenyl) 2,4 Imidazolidinedione

<

FIG-10B

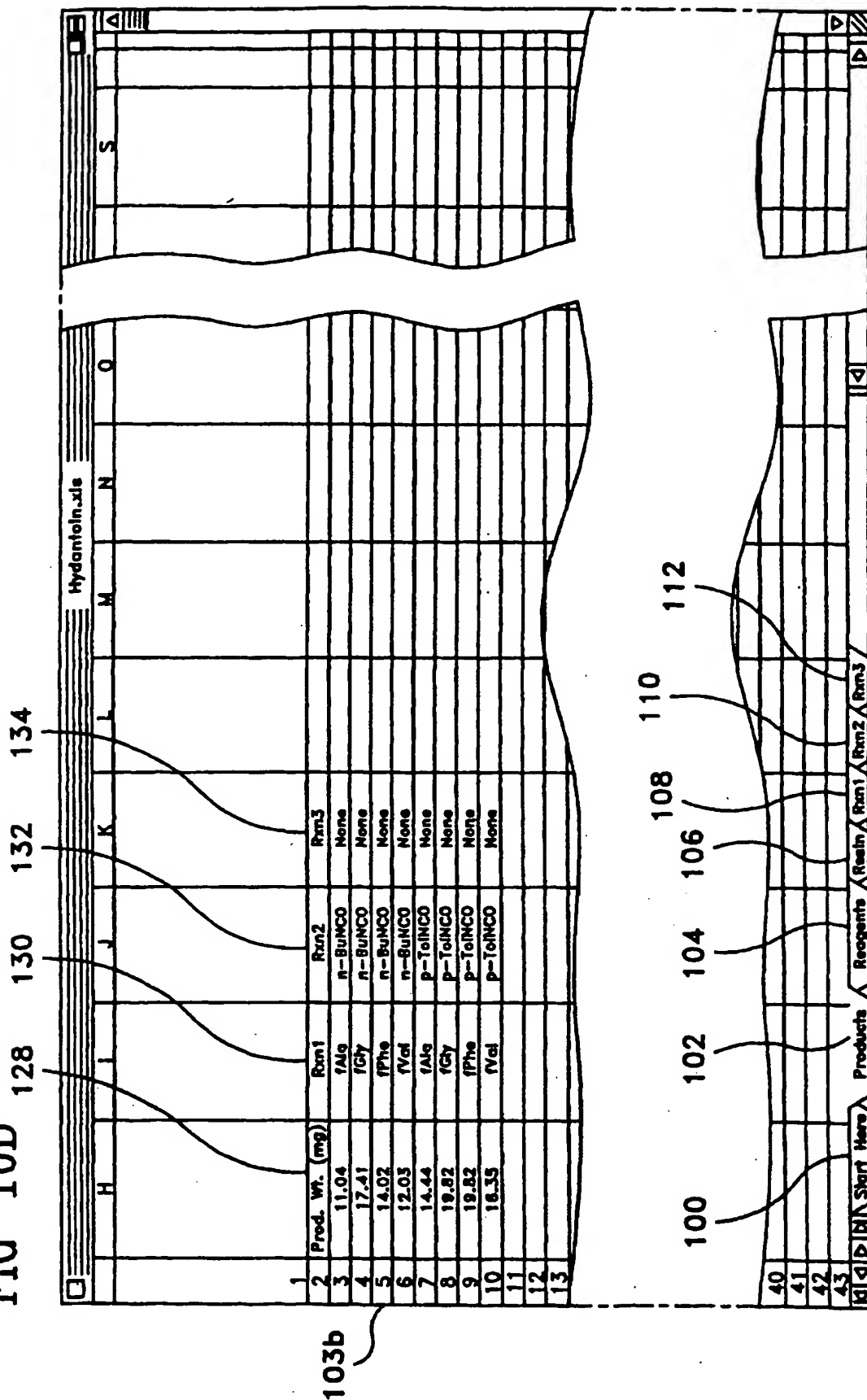


FIG-11

154

Hydantoin.xls

152

109a

100 102 104 106 108 110 112

150

164

160 162

156 158

166

168

170

172

174

176

178

180

182

184

186

188

190

192

194

196

198

200

202

204

206

208

210

212

214

216

218

220

222

224

226

228

230

232

234

236

238

240

242

244

246

248

250

252

254

256

258

260

262

264

266

268

270

272

274

276

278

280

282

284

286

288

290

292

294

296

298

300

302

304

306

308

310

312

314

316

318

320

322

324

326

328

330

332

334

336

338

340

342

344

346

348

350

352

354

356

358

360

362

364

366

368

370

372

374

376

378

380

382

384

386

388

390

392

394

396

398

400

402

404

406

408

410

412

414

416

418

420

422

424

426

428

430

432

434

436

438

440

442

444

446

448

450

452

454

456

458

460

462

464

466

468

470

472

474

476

478

480

482

484

486

488

490

492

494

496

498

500

502

504

506

508

510

512

514

516

518

520

522

524

526

528

530

532

534

536

538

540

542

544

546

548

550

552

554

556

558

560

562

564

566

568

570

572

574

576

578

580

582

584

586

588

590

592

594

596

598

600

602

604

606

608

610

612

614

616

618

620

622

624

626

628

630

632

634

636

638

640

642

644

646

648

650

652

654

656

658

660

662

664

666

668

670

672

674

676

678

680

682

684

686

688

690

692

694

696

698

700

702

704

706

708

710

712

714

716

718

720

722

724

726

728

730

732

734

736

738

740

742

744

746

748

750

752

754

756

758

760

762

764

766

768

770

772

774

776

778

780

782

784

786

788

790

792

794

796

798

800

802

804

806

808

810

812

814

816

818

820

822

824

826

828

830

832

834

836

838

840

842

844

846

848

850

852

854

856

858

860

862

864

866

868

870

872

874

876

878

880

882

884

886

888

890

892

894

896

898

900

902

904

906

908

910

912

914

916

918

920

922

924

926

928

930

932

934

936

938

940

942

944

946

948

950

952

954

956

958

960

962

964

966

968

970

972

974

976

978

980

982

984

986

988

990

992

994

996

998

1000

1002

1004

1006

1008

1010

1012

1014

1016

1018

1020

1022

1024

1026

1028

1030

1032

1034

1036

1038

1040

1042

1044

1046

1048

1050

1052

1054

1056

1058

1060

1062

1064

1066

1068

1070

1072

1074

1076

1078

1080

1082

1084

1086

1088

1090

1092

1094

1096

1098

1100

1102

1104

1106

1108

1110

1112

1114

1116

1118

1120

1122

1124

1126

1128

1130

1132

1134

1136

1138

1140

1142

1144

1146

1148

1150

1152

1154

1156

1158

1160

1162

1164

1166

1168

1170

1172

1174

1176

1178

1180

1182

1184

1186

1188

1190

1192

1194

1196

1198

1200

1202

1204

1206

1208

1210

1212

1214

1216

1218

1220

1222

1224

1226

1228

1230

1232

1234

1236

1238

1240

1242

1244

1246

1248

1250

1252

1254

1256

1258

1260

1262

1264

1266

1268

1270

1272

1274

1276

1278

1280

1282

1284

1286

1288

1290

1292

1294

1296

1298

1300

1302

1304

1306

1308

1310

1312

1314

1316

1318

1320

1322

1324

1326

1328

1330

1332

1334

1336

1338

1340

1342

1344

1346

1348

1350

1352

1354

1356

1358

1360

1362

1364

1366

1368

1370

1372

1374

1376

1378

1380

1382

1384

1386

1388

1390

1392

1394

1396

1398

1400

1402

1404

1406

1408

1410

1412

1414

1416

1418

1420

1422

1424

1426

1428

1430

1432

1434

1436

1438

1440

1442

1444

1446

1448

1450

1452

1454

1456

1458

1460

1462

1464

1466

1468

1470

1472

1474

1476

1478

1480

1482

1484

1486

1488

1490

1492

1494

1496

1498

1500

1502

1504

1506

1508

1510

1512

1514

1516

1518

1520

1522

1524

1526

1528

1530

1532

1534

1536

1538

1540

1542

1544

1546

1548

1550

1552

1554

1556

1558

1560

1562

1564

1566

1568

1570

1572

1574

1576

1578

1580

1582

1584

1586

1588

1590

1592

1594

1596

1598

1600

1602

1604

1606

1608

1610

1612

1614

1616

1618

1620

1622

1624

1626

1628

1630

1632

1634

1636

1638

1640

1642

1644

1646

1648

1650

1652

1654

1656

1658

1660

1662

1664

1666

1668

1670

1672

1674

1676

1678

1680

1682

1684

1686

1688

1690

1692

1694

1696

1698

1700

1702

1704

1706

1708

1710

1712

1714

1716

1718

1720

1722

1724

1726

1728

1730

1732

1734

1736

1738

1740

1742

1744

1746

1748

1750

1752

1754

1756

1758

1760

1762

1764

1766

1768

1770

1772

1774

1776

1778

1780

1782

1784

1786

1788

1790

1792

1794

1796

1798

1800

1802

1804

1806

1808

1810

1812

1814

1816

1818

1820

1822

1824

1826

1828

1830

1832

1834

1836

1838

1840

1842

1844

1846

1848

1850

1852

1854

1856

1858

1860

1862

1864

1866

1868

1870

1872

1874

1876

1878

1880

1882

1884

1886

1888

1890

1892

1894

1896

1898

1900

1902

1904

1906

1908

1910

1912

1914

1916

1918

1920

1922

1924

1926

1928

1930

1932

1934

1936

1938

1940

1942

1944

1946

1948

1950

1952

1954

1956

1958

1960

1962

1964

1966

1968

1970

1972

1974

1976

1978

1980

1982

1984

1986

1988

1990

1992

1994

1996

1998

2000

2002

2004

2006

2008

2010

2012

2014

2016

2018

2020

2022

2024

2026

2028

2030

2032

2034

2036

2038

2040

2042

2044

2046

2048

2050

2052

2054

2056

2058

2060

2062

2064

2066

2068

2070

2072

2074

2076

2078

2080

2082

2084

2086

2088

2090

2092

2094

2096

2098

2100

2102

2104

2106

2108

2110

2112

2114

2116

2118

2120

2122

2124

2126

2128

2130

2132

2134

2136

2138

2140

2142

2144

2146

2148

2150

2152

2154

2156

2158

2160

2162

2164

2166

2168

2170

2172

2174

2176

2178

2180

2182

2184

2186

2188

2190

2192

2194

2196

2198

2200

2202

2204

2206

2208

2210

2212

2214

2216

2218

2220

2222

2224

2226

2228

2230

2232

2234

2236

2238

2240

2242

2244

2246

2248

2250

2252

2254

2256

2258

2260

2262

2264

2266

2268

2270

2272

2274

2276

2278

2280

2282

2284

2286

2288

2290

2292

2294

2296

2298

2300

2302

2304

2306

2308

2310

2312

2314

2316

2318

2320

2322

2324

2326

2328

2330

2332

2334

2336

2338

2340

2342

2344

2346

2348

2350

2352

2354

2356

2358

2360

2362

2364

2366

2368

2370

2372

2374

2376

2378

2380

2382

2384

2386

2388

2390

2392

2394

2396

2398

2400

2402

2404

2406

2408

2410

2412

2414

2416

2418

2420

2422

2424

2426

2428

2430

2432

2434

2436

2438

2440

2442

2444

2446

2448

2450

2452

2454

2456

2458

2460

2462

2464

2466

2468

2470

2472

2474

2476

2478

2480

2482

2484

2486

2488

2490

2492

2494

2496

2498

2500

2502

2504

2506

2508

2510

2512

2514

2516

2518

2520

2522

2524

2526

2528

2530

2532

2534

2536

2538

2540

2542

2544

2546

2548

2550

2552

2554

2556

2558

2560

2562

2564

2566

2568

2570

2572

2574

2576

2578

2580

2582

2584

2586

2588

2590

2592

2594

2596

2598

2600

2602

2604

2606

2608

2610

2612

2614

2616

2618

2620

2622

2624

2626

2628

2630

2632

2634

2636

2638

2640

2642

2644

2646

2648

2650

2652

2654

2656

2658

2660

2662

2664

2666

2668

2670

2672


2674

2676

2678

2680

FIG-12

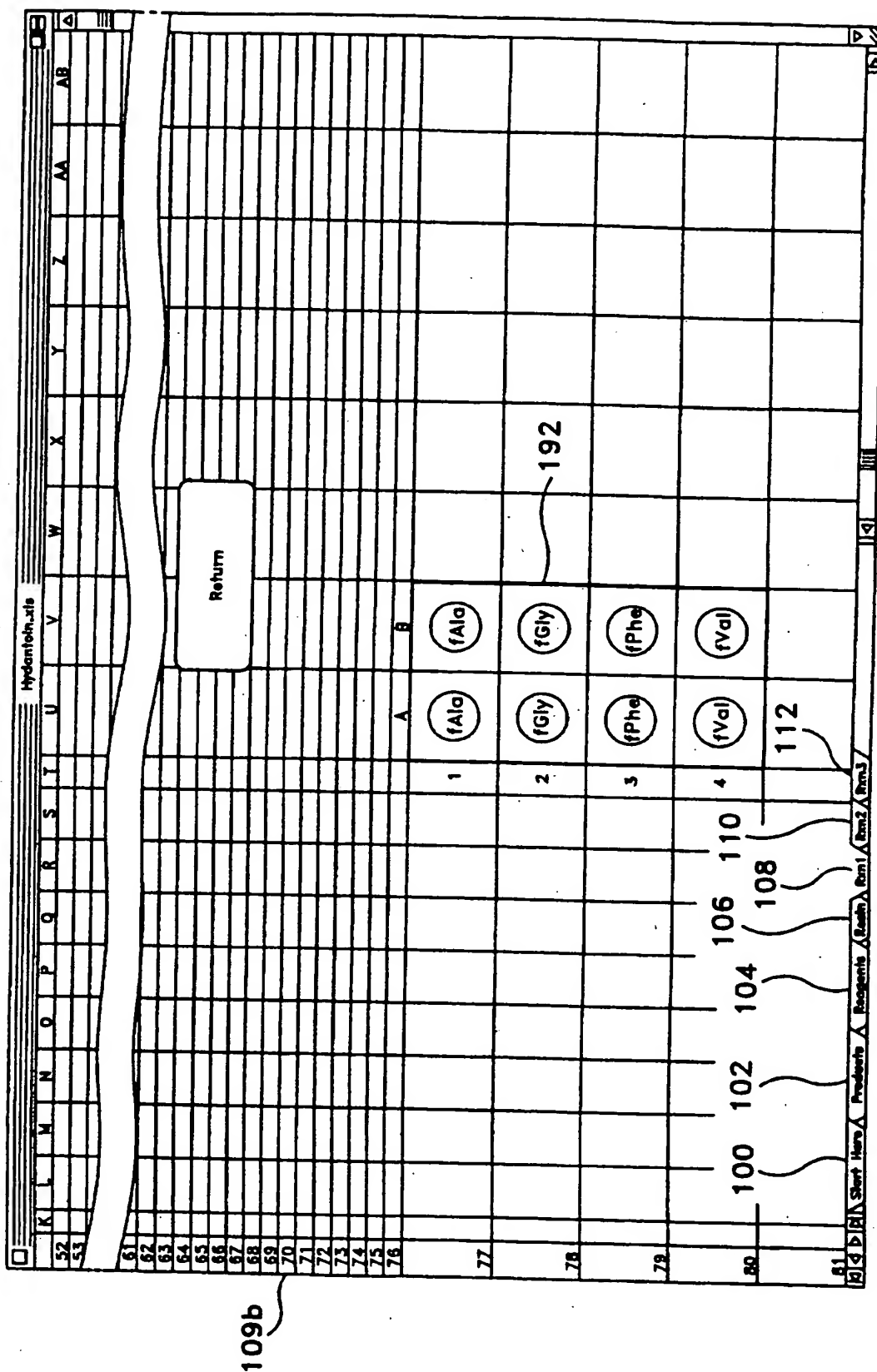



(12)

EUROPEAN PATENT APPLICATION

(88) Date of publication A3:
18.03.1998 Bulletin 1998/12

(51) Int. Cl.⁶: B01J 19/00

(43) Date of publication A2:
24.09.1997 Bulletin 1997/39

(21) Application number: 97104886.3

(22) Date of filing: 21.03.1997

(84) Designated Contracting States:
AT BE CH DE DK ES FI FR GB GR IE IT LI LU MC
NL PT SE
Designated Extension States:
AL LT LV RO SI

- **Duffield, Mark J.**
Cambridge, Massachusetts 02139-4809 (US)
- **Hogan, Eleonora M.**
Ann Arbor, Michigan 48105 (US)
- **Nickell, David Glenn**
Ann Arbor, Michigan 48103 (US)

(30) Priority: 22.03.1996 US 620265


(71) Applicant:
WARNER-LAMBERT COMPANY
Ann Arbor, Michigan 48105 (US)

(72) Inventors:
• Dewitt, Sheila Helen Hobbs
Canton, Michigan 49236 (US)

(74) Representative:
Mansmann, Ivo et al
Warner-Lambert Company,
Legal Division,
Patent Department,
c/o Gödecke AG,
Mooswaldallee 1
79090 Freiburg (DE)

(54) Information management system for automated multiple simultaneous synthesis

(57) A system (10) for developing and performing multiple simultaneous synthesis with an array of reaction vessels (32,34) includes an information management system (12), an equipment controller interface (14) and responsive equipment such as a liquid handling robot (18). Equipment may also be supplemented by other reaction peripherals (20) and pick and place robotic equipment (180) controlled by an event scheduler (22). The preferred information management system (12) includes a design module (66) accessing reagent, corporate chemical and reaction databases (68,70,72), as well as a core module (30) that includes a reaction workbook. The reaction workbook identifies the products to be produced, preferably by the product string of chemical building blocks used to form the product, the solid or solution parameters involved to generate the products, and the selection of reagents required to carry out the reactions. In addition, the reaction workbook organizes the reaction date to guide the equipment controller, preferably for movement relative to an array of multiple reaction vessels (32,34) in a reactor apparatus (28).

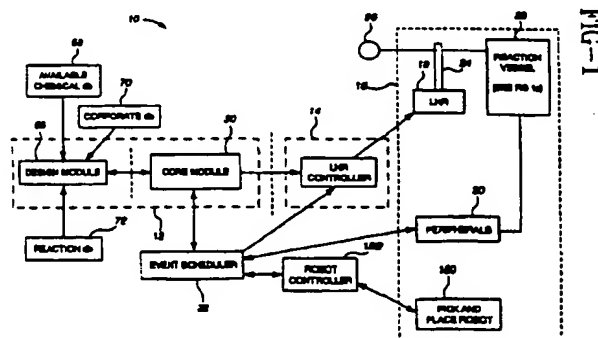

FILE-1



European Patent
Office

EUROPEAN SEARCH REPORT

Application Number
EP 97 10 4886

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.6)
X	WO 94 18226 A (UNIVERSITY OF GEORGIA RESEARCH FOUNDATION, INC.) * the whole document *	1-12, 14	B01J19/00
A	---	15-31	
A	US 5 463 564 A (DIMITRIS K. AGRAFIOTIS ET AL.) * the whole document *	1-31	
A	S.H.DEWITT ET AL.: "DIVERSOMER Technology: Solid Phase Synthesis, Automation, and Integration for the Generation of Chemical Diversity" DRUG DEVELOPMENT RESEARCH, vol. 33, no. 2, October 1994, NEW YORK, USA, pages 116-124, XP002051599 * the whole document *	1-31	
A	WO 91 17823 A (PROTOS CORPORATION) * the whole document *	1-31	TECHNICAL FIELDS SEARCHED (Int.Cl.6)
A	SUGAWARA T ET AL: "A SOLUTION-PHASE SYNTHESIS OF FRAGMENT PEPTIDE DERIVATIVES USING AN AUTOMATED SYNTHESIS APPARATUS" PROCEEDINGS OF THE SYMPOSIUM ON PEPTIDE CHEMISTRY, vol. 33, 1996, OSAKA, JP, pages 57-60, XP000617704 * the whole document *	13	B01J
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 9 January 1998	Examiner Stevnsborg, N
<p>CATEGORY OF CITED DOCUMENTS</p> <p>X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document</p> <p>T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons & : member of the same patent family, corresponding document</p>			

EPO FORM 1503 (02.92) (P/A/C/11)